

How to use the template.

1. Print pages A & B of the template and join along the line shown.
2. When joined, measure the distance for plate 23 for correct scale, this will be 23cm .
3. Cut along the outer line across the base and up the left side.
4. Turn both castor wheels toward the center (shortest distance between both).
5. Place the template on top of the existing foot plate between both castors.
6. Trim template along each numbered line, removing a section at a time until the template does not contact either castor.
7. When the plate size is determined, ensure the template fits between the lower frame rails, as the rails maybe narrower than the castors.

153


23
24
25


230

JOIN-A

CUT

CUT

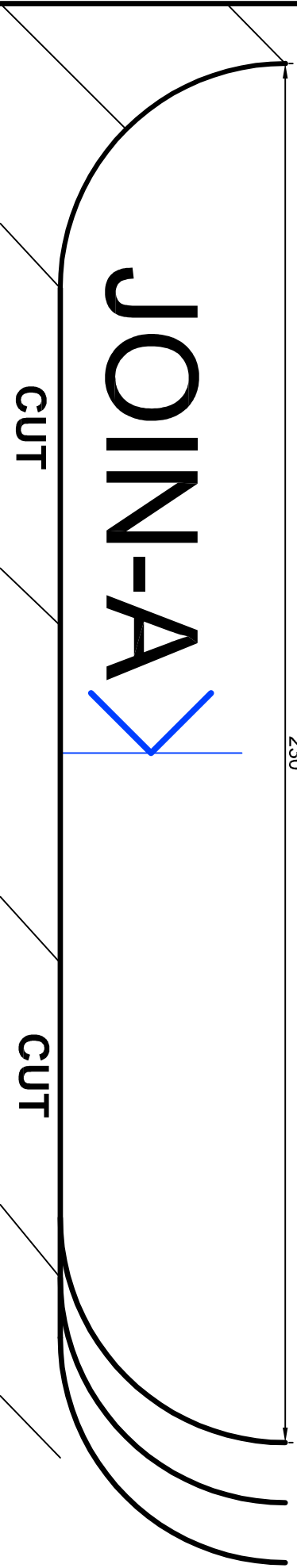

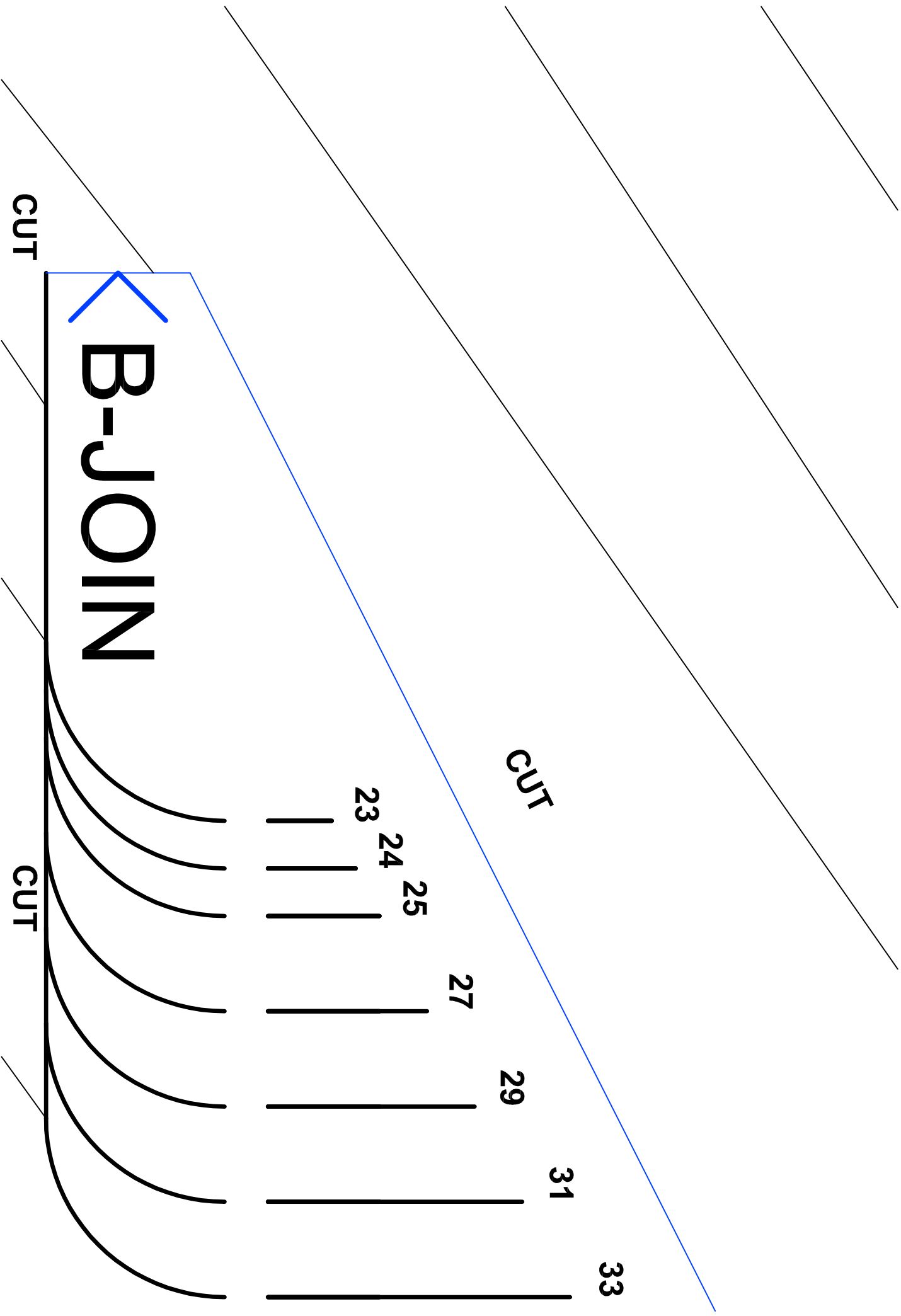

CUT

B-JOIN

CUT

CUT

23

24

25

27

29

31

33